

# J 623 2.0

## Package Inclusions

- Quick Height Adjustable (QHA) armrests
- TB3 Tilt - TRU-Balance 50 degree tilt
- Swingaway legrests with angle adjustable foot plates or centre-mount footboard
- Choice of either TRU-Comfort back or upgrade to ADI backs
- Choice of cushion
- Stealth headrest

**Further options available - refer to script form.**

With features including powerful motors that achieve speeds up to 10 km/h and improved responsiveness, the Jazzy® 623 2.0 is engineered to meet the performance needs of the most active user. The Jazzy® 623 2.0 comes standard with TRU-Balance 3 50 degree tilt and optional 25.4 cm (10") of lift elevation, while a VR2 controller offers two seating functions (as standard) and delivers smooth, nimble manoeuvrability matched with exceptional outdoor performance.

### Standard Colours



**Alternative colour options available - refer to script form.**



Attendant Control



Jazzy 623 2.0 shown in red with activated tilt



QHA Armrests



Lift




## Specifications

<b>DRIVE WHEELS</b>	Standard: 35.56 cm (14") Solid Optional: 35.56 cm (14") Pneumatic
<b>CASTOR WHEELS</b>	<b>Front:</b> 15.24 cm (6") Solid (standard) 20 cm (8") Solid (optional) <b>Rear:</b> 15.24 cm (6") Solid (standard) 20 cm (8") Solid (optional)
	N/A
<b>MAXIMUM SPEED</b> <sup>4</sup>	Up to 9.6 km/h (6 mph)
<b>GROUND CLEARANCE</b>	6.9 cm (2.7") (battery tray)
<b>TURNING RADIUS</b> <sup>5</sup>	52.1 cm (20.5")
<b>OVERALL LENGTH</b>	90.2 cm (35.5") without foot riggings
<b>BASE WIDTH</b>	62.2 cm (24.4")
<b>SEATING SIZES</b>	
<b>Synergy®/TRU-Balance® 3 Static Seat</b>	<b>Width:</b> 30 cm - 56 cm (12-22") <b>Depth:</b> 30 cm - 56 cm (12-22")
<b>TRU-Balance® 3 Power Tilt</b>	<b>Width:</b> 30.5 cm - 56 cm (12-22") <b>Depth:</b> 30.5 cm - 56 cm (12-22")
<b>SEAT-TO-FLOOR HEIGHTS</b>	
<b>Synergy®/TRU-Balance® 3 Static Seat</b> <sup>8</sup>	42 cm - 47 cm (16.5-18.5")
<b>TRU-Balance® 3 Power Tilt</b> <sup>8</sup>	44 cm - 49.5 cm (17.5-19.5")
<b>TRU-Balance® 3 Lift &amp; Tilt</b> <sup>8</sup>	44 cm - 49.5 cm (17.5-19.5")
<b>MANUAL TILT/RECLINE</b> <sup>10</sup>	Yes
<b>BATTERY SIZE</b> <sup>6</sup>	NF-22, 55 AH
<b>BATTERY WEIGHT</b> <sup>14</sup>	17.24 kg (38 lbs.)
<b>AVAILABLE ELECTRONICS</b> <sup>12</sup>	VR2 Controller
<b>BATTERY CHARGER</b>	8A Off-board
<b>MOTOR PACKAGES</b>	2-Pole Mid-Wheel 6®
<b>WEIGHT CAPACITY</b> <sup>7</sup>	136 kg (300 lbs.) <sup>9</sup>
<b>BASE WEIGHT</b>	63 kg (138.6 lbs.)



Swingaway legrests with angle adjustable foot plates or centre-mount footboard

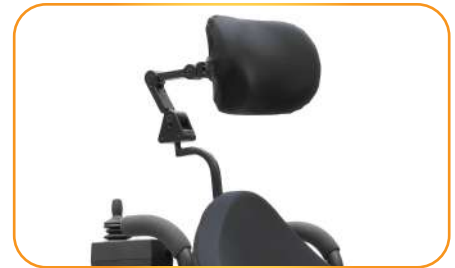

Stealth headrest

1) Standard

2) Optional

3) Requires HD Suspension

4) Speed varies with user weight, terrain type, battery charge, battery condition and tyre pressure

5) Without front riggings

6) Battery capacity may change when adding power positioning systems

7) Weight capacity is dictated by the selection of seating

8) All measurements +/- .3175 cm (1/8") with 46 cm x 46 cm (18" x 18") seat, fixed towers, 136 kg (300 lbs.) weight capacity and solid tyres (add .95 cm (3/8") with pneumatic tyres). Widths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified

9) 113 kg (250 lbs.) weight capacity with TRU-Balance® Systems

10) The Synergy® Manual Tilt weight capacity is 34 kg (75 lbs.) and maximum dimensions are 35.5 cm x 35.5 cm (14" W x 14") D. Maximum dimensions of the Manual Recline vary by power base and seating system

11) Requires Group 34 batteries

12) All Quantum Rehab® units can be ordered with specialty controls

13) Required with high-speed motors

14) Battery weight may vary +/- 1.4 kg (3 lbs.) based on manufacturer

15) Call for availability